


AlphaLoc

Overdenture Attachments

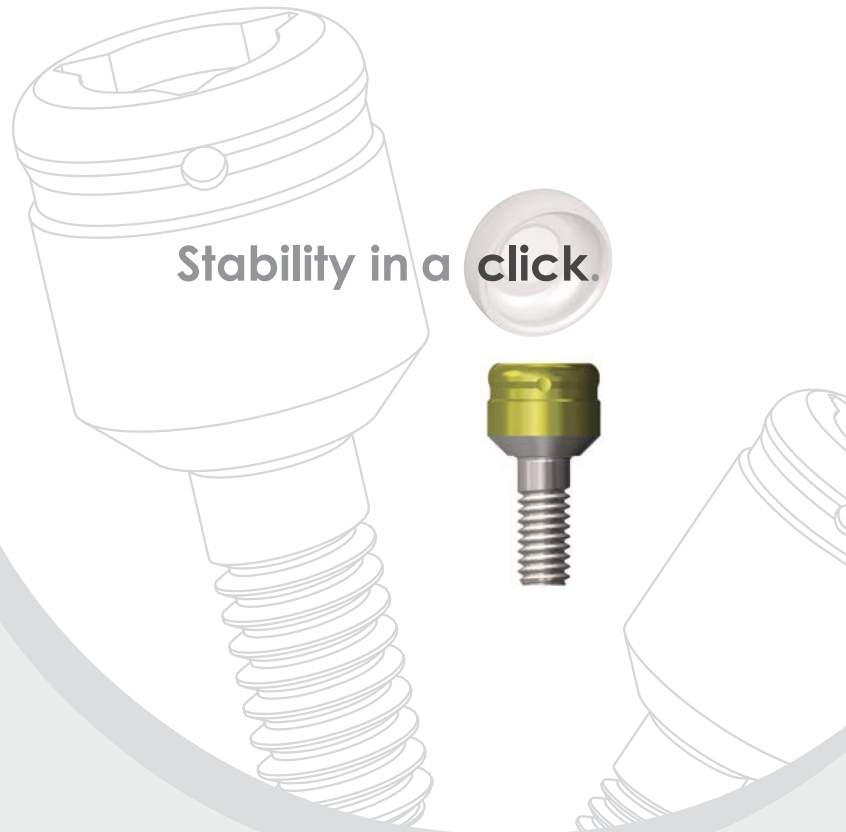

 AlphaBio^{TEC}
Implantology

Stability in a click.

The AlphaLoC Locator Overdenture Attachments is an effective, simple and easy solution for overdenture restoration. Alpha-Bio Tec's AlphaLoC Locator system has been especially designed for optimal compatibility with Alpha-Bio Tec's internal hex implant systems for maximal restoration precision.

The AlphaLoC Locator is simple for the patient. The self locating attachments enable easy and correct overdenture seating therefore minimizing attachment wear caused by improper overdenture seating.

The AlphaLoC Locator is simple for the dentist. The system offers a variety of options for gingival height, retention and angle correction for a wide range of cases. Better overdenture positioning means less wear of the nylon liners, therefore fewer patient recall visits.

AlphaLoC Locator advantages:

Low vertical height

- Enables treatment of a wide range of cases.
- Abutment heights vary from as low as 0.3 mm (2.5 mm including attachment) – ideal for use in tight interocclusal spaces.

Simple and easy to use

- Overdenture attachment simply guides itself to the correct place on the abutment
- Accurate overdenture placement

Reliable, long term solution

- Dual retention: a combination of inside and outside retention ensures a durable and long lasting solution
- Greater retention surface area compared to internal retention attachments
- The design of the pivoting AlphaLoC Locator Male allows a resilient connection for the overdenture while maintaining optimal retention

Versatile treatment options

- A variety of abutment heights for maximal suitability to gingival height
- A variety of angles:
 - Up to 20° using AlphaLoC Locator Abutments
 - Up to 50° using AlphaLoC Locator UniCover
- A variety of retention levels depending on the chosen replacement male:
 - Regular range - three retention levels: 5.0 lbs. (clear), 3.0 lb. (pink), 1.5 lb. (blue)
 - Extended range - three retention levels: 4.0 lbs. (green), 2.0 lb. (orange), 1.0 lb. (red)

AlphaLoC Locator Abutment

Titanium alloy with TiN coating



Height	0.3 mm	1 mm	2 mm	3 mm	4 mm	5 mm	6 mm
Ref No.	4801	4802	4803	4804	4805	4806	4807

AlphaLoC Locator UniCover

Titanium alloy with TiN coating



Height	1 mm	2 mm
Ref No.	4833	4834

For Alpha Universe Multi Unit Abutment System

Catalog Numbers

AlphaLoC Locator Male Processing Package

Titanium Cap, Block Out Spacer, Black Low Density Polyethylene Male, Nylon Replacement Males (clear, pink, blue)



Content	1 unit
Ref No.	4811

AlphaLoC Locator Extended Range Male Processing Package

Titanium Cap, Block Out Spacer, Black Low Density Polyethylene Male, Nylon Extended Range Replacement Males (green, orange, red)



Content	1 unit
Ref No.	4812

AlphaLoC Locator Kit 0° - 10°

AlphaLoC Locator Abutment
AlphaLoC Locator Male Processing Package



Height	0.3 mm	1 mm	2 mm	3 mm	4 mm	5 mm	6 mm
Ref No.	4845	4846	4847	4848	4849	4850	4851

Select required kit according to abutment height

AlphaLoC Locator Kit 10° - 20°

AlphaLoC Locator Abutment
AlphaLoC Locator Extended Range Male Processing Package



Height	0.3 mm	1 mm	2 mm	3 mm	4 mm	5 mm	6 mm
Ref No.	4852	4853	4854	4855	4856	4857	4858

Select required kit according to abutment height

AlphaLoC Locator Replacement Male



Content	4 units	4 units	4 units
Ref No.	4836	4837	4838

AlphaLoC Locator Extended Range Replacement Male



Content	4 units	4 units	4 units
Ref No.	4839	4840	4841

AlphaLoC Locator Black Processing Replacement Male



Content	4 units
Ref No.	4842

AlphaLoC White Block Out Spacer



Content	4 units
Ref No.	4843

AlphaLoC Locator Impression Coping



Content	4 units
Ref No.	4830

AlphaLoC Locator Female Analog (4 mmd)



Content	4 units
Ref No.	4835

AlphaLoC Locator Parallel Post



Content	4 units
Ref No.	4844



Content	1 unit
Ref No.	4825

AlphaLoC Locator Core Tool



Content	1 unit	1 unit	1 unit
Ref No.	4808	4809	4810

AlphaLoC Locator Male Removal Tool (Tip End Only)



AlphaLoC Locator Abutment Driver (Gold End Only)



AlphaLoC Wrench Locator Driver (15 mm)



Content	1 unit
Ref No.	4831


AlphaLoC Locator Motor Mount Driver (23 mm)



Content	1 unit
Ref No.	4832



Stability in a click.

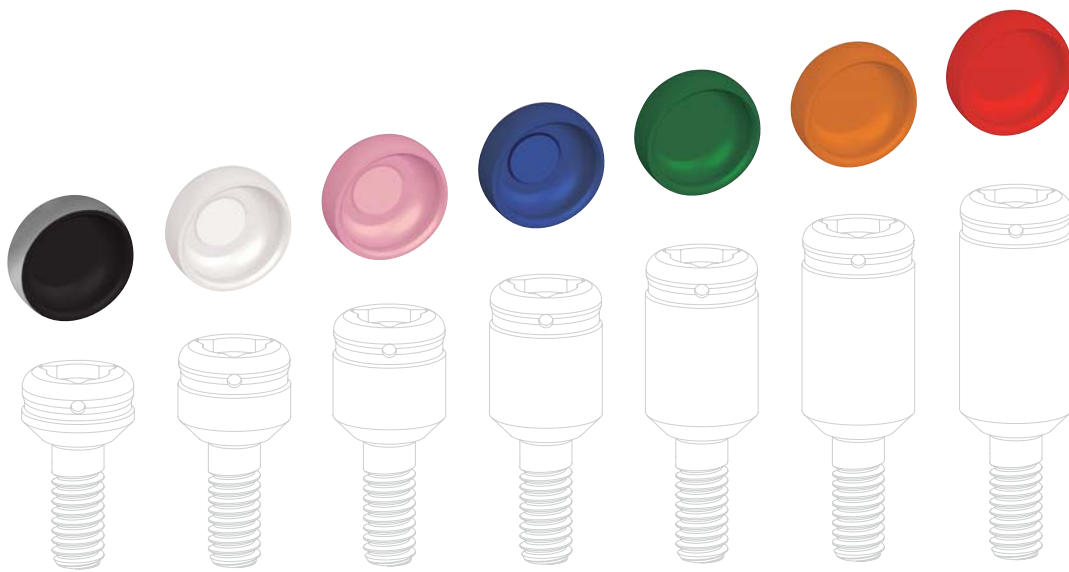

 **ALPHA BIO**^{TEC}
 Implantology

www.alpha-bio.net

Alpha-Bio Tec's products are cleared for marketing in the USA and are CE-marked in accordance with the Council Directive 93/42/EEC and Amendment 2007/47/EC. Alpha-Bio Tec complies with ISO 13485:2003 and the Canadian Medical Devices Conformity Assessment System (CMDCAS). Product availability may vary between countries.

Alpha-Bio Tec Ltd.

7 Hatnufa St., Kiryat Arye
 P.O.B. 3936, Petach Tikva
 49510, Israel

Israel

T+972-3-9291000
 F+972-3-9235055
sales@alpha-bio.net

International

T+972-3-9291055
 F+972-3-9291010
export@alpha-bio.net

MEDES LIMITED

5 Beaumont Gate, Shenley Hill,
 Radlett, Herts WD7 7AR,
 England
 Tel / Fax: +44 1923859810